

## Linking made easy

Via a QR code



Comfortable  
and easy

Connect the freeze dryer  
with your smartphone  
via QR code



Before using it for the first time, the freeze dryer and a smartphone must be linked. For this purpose, a unique QR code is generated in the freeze dryer during the initialisation. This QR code has to be scanned with your smartphone. This ensures that the process data can only be viewed on authorised smartphones. The link can be undone any time via the freeze dryer.



## Stay connected



LyoApp –  
connecting your  
freeze dryer with  
your smartphone

**Martin Christ**  
**Gefriertrocknungsanlagen GmbH**  
An der Unteren Söse 50  
37520 Osterode • Germany  
Tel. +49 (0) 55 22 50 07-0  
Fax +49 (0) 55 22 50 07-12  
info@martinchrist.de  
www.martinchrist.de


## For greater independence

### LyoApp

The LyoApp is an innovative tool for linking your pilot freeze dryer to your smartphone. The app is particularly interesting for customers who would like to stay continuously informed about the freeze drying process in any type of situation.

#### The LyoApp offers the following features:

- Overview of the current process data
- Name of the current process section
- Zoomable visualisation of the process graphics
- Important status information provided as push notifications
- Messages and alerts
- Configuration of shown process data
- Linking with several freeze dryers

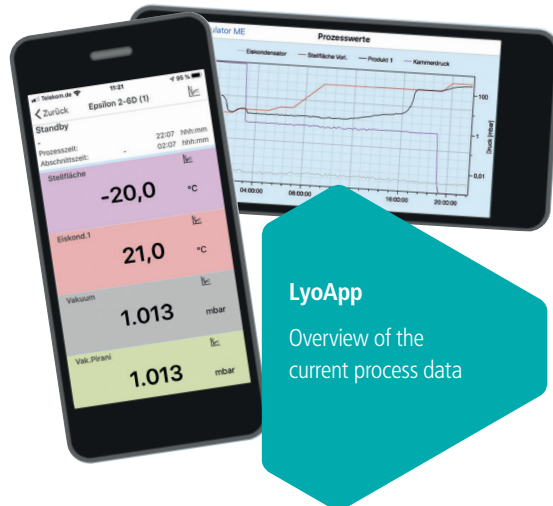

## Uncompromised security

### Secure encryption



The process data are transferred to the LyoApp of the connected smartphones via the Christ Internet Service.

The communication is entirely wireless. The freeze dryer information is encrypted before the transfer so that it can only be accessed by the connected smartphones. To use this feature, the freeze dryer must have access to the internet.

## Benefits

### Wide range of smart features

Users can benefit from the LyoApp and its extensive features in numerous ways.

#### Advantages at a glance:

- Transparent overview of the freeze drying process
- Messages and alerts which enable the user to react whenever necessary
- Can be used anytime and anywhere – from the office next door or from the other end of the world:
  - via intranet
  - via internet
- High level of security thanks to data encryption
- No process data storage in the Christ Internet Services
- Process data are constantly available, including after interruptions
- Available for current operating systems of Android and iPhone smartphones

Currently available  
for new ordered  
pilot freeze dryers

Epsilon 2-6D LSCplus and  
Epsilon 2-10D LSCplus